blue●dot®

B-FT102 BT

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EN : BLUEDOT SMART INFRARED THERMOMETER

Thank you for your purchase. This product is designed specifically to measure/monitor body temperature for homecare and medical use. Please read all instructions carefully before using, and retain this manual for future reference.

Table of Contents

1. Key of symbols	02
2. Product Features	05
3. Product Information	05
1.Product appearance.....	05
2.Lattice screen display instructions.....	05
3. Error and Troubleshooting.....	06
4. Product Specifications	07
5. User Instructions	07
1. How to connect the app.....	07
2. Use method.....	08
3. Fever alarm.....	09
4. Memory and mute function.....	10
5. Unit conversion.....	10
6. Turn off.....	10
7. Battery installation.....	10
8. Storage & Maintenance.....	11
6. Included Contents	11
7. Classification	12
8. Additional Information	12
9. Disposal	13
10. Warranty/service	13

1. Key of symbols

Symbol	Description
 WARNING	WARNING indicates a hazardous situation which, if not avoided, could result in death or serious injury
 CAUTION	CAUTION indicates a hazardous situation which, if not avoided, may result in minor or moderate injury
 NOTICE	NOTICE addresses practices not related to personal injuries, such as product and/or property damage
	Type BF
	Do not dispose of this product as unsorted municipal waste. Collection of such waste separately for special treatment is necessary
 0123	This item is compliant with the Medical Device Directive 93/42/EEC of June 14, 1993, a directive of the European Economic community
	Refer to instruction manual/booklet
IP22	Waterproof degree
	Manufacturer
	Serial number
	Keep Dry
	Fragile
	Upside
	Maximum stacking layer

 **WARNING**

- If you are concerned about whether you are having a fever or not, please consult a physician. This device is no substitute for medical treatment.
 - Please dispose the batteries and device in accordance with local regulations.
 - The device is not suitable for use in environment of flammable anesthetic mixtures.
 - The device requires no calibration.
 - The device is not repairable and contains no repairable parts. No modification of this equipment is allowed.
 - Please do not dismantle, re-assemble the device, and do not fall the device heavily, and avoid strong impact and shock.
 - The manufacturer need no preventive inspections by others.
 - The thermometer is used only for personal monitoring, it cannot replace a doctor.
 - Please do not use the thermometer for other purpose.
 - When using on children, please observe the general safety precautions.
 - Please take out the batteries if the device is unused for a long time(more than 3 months).
- If the device has been damaged or there is abnormal, please stop using.
- Incorrect usage of battery may cause damage.
 - Please do not re-charge or dismantle the battery or put the battery near or into the fire, to avoid the battery explosion.
 - Please do not operate the device in environment out of the operating temperature range.
 - Please avoid biting, throwing, folding, pulling or any strong impact on the device, or it may cause failure while measuring.

- Keep the device away from kids.

 **CAUTION**


- Please do not dip the device into water or other liquids.
- Please follow the maintenance and storage section in the manual when cleaning or disinfecting.
- Keep the device away from direct sunlight and store it in dry and dust-free environment, ideally within storage temperature range of -20°C to $+55^{\circ}\text{C}$ (-4°F to $+131^{\circ}\text{F}$).
- Please do not use if the sensor or product is broken. Do not repair it by oneself if the device is broken, please contact your nearest customer service center.
- Please make sure the sensor mirror is clean before measurement, or please clean it with alcohol cotton before use.
- Please avoid dropping or strong impact to protect the precision parts.
- Please do not twist the device or sensor.
- Please wait for at least 30 seconds after measuring each time to make sure accurate result.
- Please wait for a few minutes until the sensor return to its normal temperature, after cleaning the sensor with alcohol.
- It's recommended to wait for a few minutes in the same environment until baby's body temperature becomes stable if drastic temperature changes are observed, as baby are slower at adapting to temperature changes in the environment.
- It is recommended to take 3 times of measurement, and choose the highest for reference for babies at the age of 3 months to 3 years, or for first-time users of the device.

2. Product Features


- FDA cleared for home and professional use
- Instant temperature reading
- 32-memory function to track fever
- Fever alert with different backlight color
- Modern and sleek design with large backlit LED display for clear visibility in darkness
- Fast °C/°F temperature reading

3. Product Information

1. Product appearance

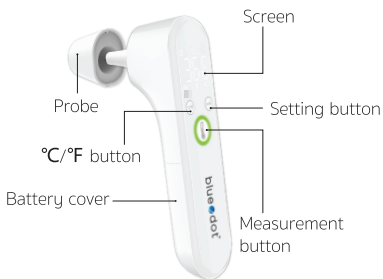2. Lattice screen display instructions

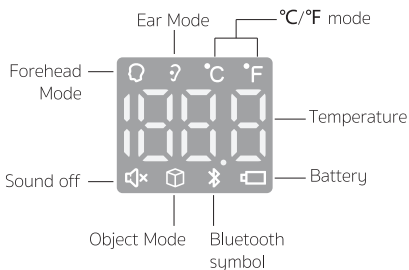3. Error and Troubleshooting

Symptom	Possible Cause	Description & Solution
	Low Battery	Replace batteries.
H	Results The measured human body temperature was high.	Check whether it is aimed at human body parts.
L	Results The human body temperature was low.	Check whether it is aimed at human body parts.
EH	The environment is too high.	Check whether the ambient temperature is within the specification, and move the product to the ambient temperature required by the specification for half an hour before measuring.
EL	The environment is too low.	Check whether the ambient temperature is within the specification, and move the product to the ambient temperature required by the specification for half an hour before measuring.
Er A	The temperature is unstable (the default shutdown does not show)	When the temperature changes, please place the product in the current environment for half an hour before use.
Err	system malfunction	Please place the product and try again. If it is still abnormal, please contact your local dealer.
No display.	The battery is exhausted.	Please replace the battery.

4. Product Specifications

Product type	IR Thermometer
Model No.	B-FT102 BT
Measurement range	32.0°C ~ 43.0°C (89.6°F ~ 109.4°F)
Accuracy	35.0°C ~ 42.0°C, ±0.2°C (±0.4°F) Within other ranges, 0.3°C (0.5 °F)
Resolution	0.1°C/0.1°F
Unit	°C/°F
Measurement time	1 second
Power supply	2xAAA(1.5V), DC 3.0V
Working environment	Temperature: 10°C ~ 40°C Relative humidity: ≤85% RH Atmospheric pressure: 70.0kPa ~ 106.0kPa
Storage environment	Temperature: -20°C ~ +55°C Relative humidity: ≤93% RH Atmospheric pressure: 70.0kPa ~ 106.0kPa
Product size	140(L)×36(W)×42(H) mm.
Weight	66 g (excluding battery)

5. User Instructions

1. How to connect the app

1. Download the free “Smart bluedot” application from the App Store (IOS) or from Google Play (Android).
2. Open the application “Smart bluedot”.
3. The system will prompt for registration. for you to proceed with registration (If you have already registered. Please skip to step 5).

4. Complete the information as per the steps in the application.
***Required to specify height, gender, day-month-year (year) of birth**
5. When the data is saved successfully, log-in into the system, select the

6. Once the device is paired successfully. You can immediately measure body temperature. when the measurement is completed The measured value will be displayed on the application. Then press record the measured values.

2. Use method

Temperature measurement in forehead mode:

- Press the measurement key for a short time.

and then turn it on for full display before standby (in standby, the icons "°C or "°F" and "🔋" are always on without flashing)

- Please keep the probe 3-5 cm from the center of the user's forehead, then press the measurement button. A beep indicates the measurement has been finished.

- The device will automatically turn off without any operation or you can press the measurement button for 3 seconds to turn off the power.

Note:

- Please store the device in a proper environment. Or place the device in ambient temperature within the range of 10°C ~ 40 °C for at least 30 minutes in advance, to avoid inaccurate results due to extreme cold or hot.
- Take a rest for at least 30 minutes before measurement, after strenuous exercise.
- Make sure there is no hair, sweat, hat, or makeup on your forehead during the temperature measurement.

- It is recommended for users to stay in the same room for at least 20 minutes to stabilize the body temperature before measurement for the accurate results.
- It is recommended to take multiple measurements to get the accurate result to avoid human error.

Ear temperature measurement:

- Take off the probe cover (it is automatically switched to ear mode), and fit the probe snugly into the ear canal. and then turn it on for full display before standby (in standby, the icons "°C or °F" and "? " are always on without flashing)
- Press and release the Measure/Mute-unmute button in 1 second, the beep is heard, you can now read the value.
- The device will automatically turn off without any operation or you can press the measurement button for 3 seconds to turn off the power.

Object temperature measurement:

- When the machine is turned off, press and hold the "memory/mute key" and the "measure key" for 3 seconds at the same time, let go and enter the object temperature mode.

Temperature range of object: 0-93.2°C(0-199.9°F)

- At this time, the icon lights up, and you can start to click the measure key to measure the object temperature.
- The device will automatically turn off without any operation or you can press the measurement button for 3 seconds to turn off the power.

3. Fever alarm :

When the measurement is higher than the alarm value, the device will emit a continuous beep sound and the button will turn to different backlight color to remind the user to pay attention.

Temperature range	Button backlight color	Beep mode
<37.5°C(99.5°F)	Green	One short beep
37.5°C~38.5°C (99.5°F~101.3°F)	Yellow	2 times of double short beep "Beep-beep...Beep-beep...Beep-beep..."
>38.5°C(>101.3°F)	Red	3 times of double short beep "Beepbeepbeep ... Beepbeepbeep... Beepbeepbeep..."

4. Memory and mute function :

- Short Press the the “memory/mute key” button to check memory instand by status.
- Hold the “memory/mute key” button to switch to mode in mode standbystatus.
- Take out the batteries to cut off power supply and the memory will be clean automatically.

5. Unit conversion:

Short Press the the“Mode” button to switchbetween °C and °F, and the default unit is °C.

6. Turn off

- Manual power off: Inthe standbystatus, press the measurement button for three seconds.
- Automatic shut-down : 15 seconds automatically power off without any operation.

7. Battery installation

Please replace the batteries when the low battery displays on the screen, prepare 2 AAA batteries and install as follows:

- Push the battery cover at the bottom of the device along the arrow direction.
- Take out the old batteries and install the new ones correctly according to the positive “+” and negative “-” symbol in the bottom. Then close the battery cover.
- Please take out the batteries to prevent battery leakage, if used for more than 3 months.

8. Storage & Maintenance

- Wipe with a dry, soft cloth for daily cleaning. Wipe with water or neutral detergent first if contamination is severe, and then wipe it with a dry cloth.
- Please take out the batteries to prevent battery leakage, if the device has not been used for more than 3 months.
- It is recommended to clean with 75% alcohol before and after each use. This product does not require disinfection.
- Keep the device away from direct sunlight and store in dry and dust-free environment, relative humidity should be lower than 93%, temperature should be +5 ~ +40°C.
- Please do not use petrol, volatile oil, thinner, alcohol and other chemical solvents to wipe the product. If necessary, wipe with a damp cloth moistened with water, but avoid the water entering.

6. Included in delivery

1.5 V AAA batteries	2 pcs.
Instruction manual (Eng/Thai)	1 pc.
IR forehead thermometer	1 pc.

7. Classification

Internally powered equipment;

Type BF applied part: Probe;

IP22;

Not category AP / APG equipment;

Mode of operation: continuous operation.

8. Additional Information

Non-toxic, no stimulation, no sensitization adverse reaction for the user during the test.

NOTICE

- The B-FT102 BT IR forehead thermometer meets YY0505 standard electromagnetic compatibility requirements.
- Users should install and use according to the electromagnetic compatibility information provided by the manual.

NOTICE

- Portable and mobile RF communication equipment may affect the performance of device. Avoid getting close to strong electromagnetic interference when using it, such as mobile phones, microwave ovens and etc.
- The B-FT102 BT IR forehead thermometer should not be used close to or stacked with other equipment. If you have to please check if it can operate normally before use.

9. Disposal

Empty, completely flat batteries must be disposed of through specially designated collection boxes, recycling points or electronics retailers. You are legally required to dispose of the batteries.

The codes below are printed on batteries containing harmful substances:

- Pb = Battery contains lead
- Cd = Battery contains cadmium
- Hg = Battery contains mercury

For environmental reasons, do not dispose of the device in the household waste at the end of its service life.

Dispose of the device at a suitable local collection or recycling point in your country.

Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.

Dispose of packaging in an environmentally friendly manner

10. Warranty/service

LIFEBOX ASIA Co.,Ltd warrantys this product to be materially free from defects and workmanship to the original purchaser for a period of 2 years from the original purchase date.

LIFEBOX ASIA Co.,Ltd responsibility with regard to this limited warranty shall be limited solely to repair or replacement at its option, of any product which fails during normal consumer use.

This warranty does not extend to damage or failure which results from misuse, neglect, accident, alteration, abuse, improper installation or maintenance. At any time during the 1 years period following the purchase, the product fails due to defects in material or workmanship, return the defective product with proof of purchase.

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